

## Carcinoembryonic Antigen (CEA) / CD66 Antibody

Mouse Monoclonal Antibody [Clone C66/1009]


Catalog No	Format	Size
1048-MSM5-P0	Purified Ab with BSA and Azide	200ug/ml
1048-MSM5-P1	Purified Ab with BSA and Azide	200ug/ml
1048-MSM5-P1ABX	Purified Ab WITHOUT BSA and Azide	1.0mg/ml

Applications	Tested Dillution
Flow Cytometry (Flow)	1-2ug/million cells
Immunofluorescence (IF)	1-3ug/ml
Immunohistochemistry (IHC)	1-2ug/ml


Product Details	
Clone	C66/1009
Gene Name	CEACAM5
Immunogen	Recombinant full-length human CEA protein
Host	Mouse
Clonality	Monoclonal
Isotype / Light Chain	IgG2a
Mol. Weight of Antigen	80-200kDa
Cellular Localization	Apical cell membrane, Cell membrane, Cell surface
Species Reactivity	Human
Positive Control	MCF7 or 293T cells. Human colon carcinoma.

*\*Optimal dilution for a specific application should be determined.*

### Product Images for Carcinoembryonic Antigen (CEA) / CD66 Antibody

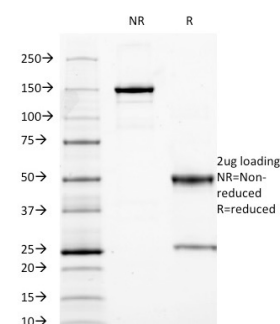

SDS-PAGE Analysis of Purified CEA Mouse Monoclonal Antibody (C66/1009). Confirmation of Purity and Integrity of Antibody.

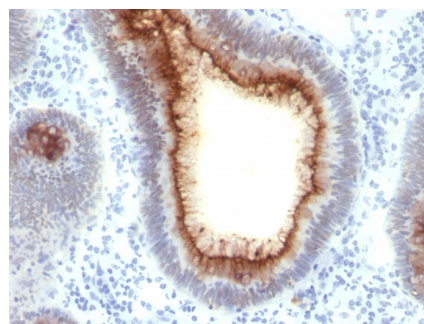

Formalin-fixed, paraffin-embedded human colon carcinoma stained with CEA Mouse Monoclonal Antibody (C66/1009).

### Specificity & Comments

This antibody recognizes proteins of 80-200kDa, identified as different members of CEA family. CEA is synthesized during development in the fetal gut and is re-expressed in increased amounts in intestinal carcinomas and several other tumors. This MAb does not react with nonspecific cross-reacting antigen (NCA) and with human polymorphonuclear leucocytes. It shows no reaction with a variety of normal tissues and is suitable for staining of formalin/paraffin tissues. CEA is not found in benign glands, stroma, or malignant prostatic cells. Antibody to CEA is useful in detecting early foci of gastric carcinoma and in distinguishing pulmonary adenocarcinomas (60-70% are CEA+) from pleural mesotheliomas (rarely or weakly CEA+). Anti-CEA positivity is seen in adenocarcinomas from the lung, colon, stomach, esophagus, pancreas, gallbladder, urachus, salivary gland, ovary, and endocervix.

---

### Research Areas

Hematopoietic Stem Cells

---

### Known Applications & Suggested Dilutions

Flow Cytometry (0.5-1ug/million cells) | Immunofluorescence (1-2ug/ml) | Immunohistochemistry (Formalin-fixed) (1-2ug/ml for 30 minutes at RT), (Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95 &degC followed by cooling at RT for 20 minutes), Optimal dilution for a specific application should be determined.

---

### Limitations and Warranty

This antibody is available for research use only and is not approved for use in diagnosis.

There are no warranties, expressed or implied, which extend beyond this description. Company is not liable for any personal injury or economic loss resulting from this product.

---

### Supplied As

200ug/ml of Ab purified from Bioreactor Concentrate by Protein A/G. Prepared in 10mM PBS with 0.05% BSA & 0.05% azide. Also available WITHOUT BSA & azide at 1.0mg/ml.

---

### Storage and Stability

Antibody with azide - store at 2 to 8 °C. Antibody without azide - store at -20 to -80 °C. Antibody is stable for 24 months. Non-hazardous. No MSDS required.

---