

# Ki-67 (Proliferating Cell Marker) Antibody

Mouse Monoclonal Antibody [Clone MKI67/2461]


Catalog No	Format	Size
4288-MSM1-P0	Purified Ab with BSA and Azide at 200ug/ml	20 ug
4288-MSM1-P1	Purified Ab with BSA and Azide at 200ug/ml	100 ug
4288-MSM1-P1ABX	Purified Ab WITHOUT BSA and Azide at 1.0mg/ml	100 ug


Applications	Tested Dillution	Note
Flow Cytometry (Flow)	1-2ug/million cells	
Immunofluorescence (IF)	1-3ug/ml	
Immunohistochemistry (IHC)	1-2ug/ml	30 min at RT. Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95°C followed by cooling at RT for 20 minutes
Western Blot (WB)	2-4ug/ml	


MKI67/2461	
MKI67	
Recombinant fragment of human Ki67 protein (around aa 2293-2478) (exact sequence is proprietary)	
Mouse	
Monoclonal	
IgG1 / Kappa	
345kDa and 395kDa	
Chromosome, Nucleolus, Nucleus	
Human	
Any actively proliferating cells. Human skin, tonsil or lymph node.	

<sup>\*</sup>Optimal dilution for a specific application should be determined.


# Product Images for Ki-67 (Proliferating Cell Marker) Antibody

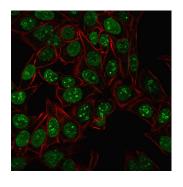


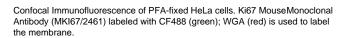

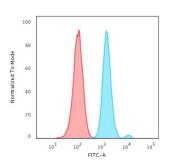



Flow Cytometric Analysis of HeLa cells using Ki67 Mouse Monoclonal Antibody (MKI67/2461) followed by goat anti-mouse IgG-CF488 (blue); isotype control (red)

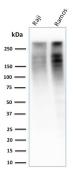


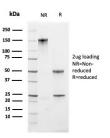

Western Blot Analysis of Raji and Ramos cell lysates using Ki67 Mouse Monoclonal Antibody (MKI67/2461).

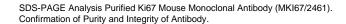


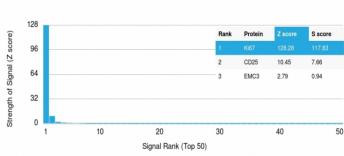

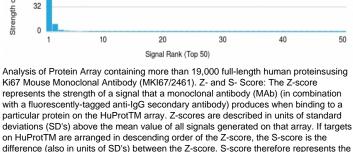

relative target specificity of a MAb to its intended target. A MAb is considered to specific to its intended target, if the MAb has an S-score of at least 2.5. For example, if a MAb binds to protein X with a Z-score of 43 and to protein Y with a Z-score of 14.

then the S-score for the binding of that MAb to protein X is equal to 29.

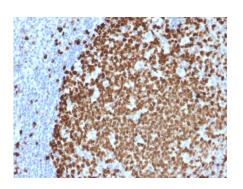

Formalin-fixed, paraffin-embedded human tonsil stained with Ki67 Mouse Monoclonal Antibody (MKI67/2461).

### **Specificity & Comments**

Ki-67 antigen is a nuclear, non-histone protein that is present in all stages of the cell cycle except G0. This characteristic makes Ki-67 an excellent marker for proliferating cells and is commonly used as one of the prognostic factors in cancer studies. A correlation has been demonstrated between Ki-67 index and the histo-pathological grade of neoplasms. Assessment of Ki-67 expression in renal and ureter tumors shows a correlation between tumor proliferation and disease progression, thus making it possible to differentiate high-risk patients. Ki-67 expression may also prove to be important for distinguishing between malignant and benign peripheral nerve sheath tumors. Ki-67 labeling index has been shown to be a prognostic marker in a number of neoplasms including grade II astrocytoma, oligodendroglioma, colon carcinoma, and breast carcinoma. In general, Ki-67 is a good marker of proliferating cell populations.

#### Supplied As

200ug/ml of Ab purified from Bioreactor Concentrate by Protein A/G. Prepared in 10mM PBS with 0.05% BSA & 0.05% azide. Also available WITHOUT BSA & azide at 1.0mg/ml.

## Storage and Stability

Antibody with azide - store at 2 to 8 °C. Antibody without azide - store at -20 to -80 °C. Antibody is stable for 24 months. Non-hazardous. No MSDS required.

#### Research Areas

Neuroscience, Nuclear Marker

#### **Limitations and Warranty**

This antibody is available for research use only and is not approved for use in diagnosis. There are no warranties, expressed or implied, which extend beyond this description. Company is not liable for any personal injury or economic loss resulting from this product.

