

Cytokeratin, Multi (Epithelial Marker) Antibody

Mouse Monoclonal Antibody [Clone C11]

Catalog No	Format	Size
MSM1-393-P0	Purified Ab with BSA and Azide	200ug/ml
MSM1-393-P1	Purified Ab with BSA and Azide	200ug/ml
MSM1-393-P1ABX	Purified Ab WITHOUT BSA and Azide	1.0mg/ml


Applications	Tested Dillution
Flow Cytometry (Flow)	1-2ug/million cells
Immunofluorescence (IF)	1-3ug/ml
Immunohistochemistry (IHC)	1-2ug/ml
Western Blot (WB)	2-4ug/ml

Product Details

C11	
KRT4	
Keratin-enriched preparation from cultured human A431	
Mouse	
Monoclonal	
IgG1 / Kappa	
Multiple	
Cow, Frog, Goat, Guinea Pig, Human, Marmoset, Mouse, Pig, Rat	
HeLa or A431 cells. Skin or colon carcinoma.	

*Optimal dilution for a specific application should be determined.

Product Images for Cytokeratin, Multi (Epithelial Marker) Antibody

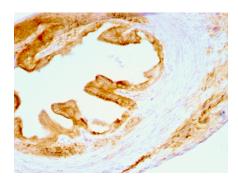


 $250 \rightarrow$ $150 \rightarrow$ $100 \rightarrow$ $75 \rightarrow$ $37 \rightarrow$ 22g loading NR=Non-reduced R=reduced R=reduced 15 \rightarrow $10 \rightarrow$

NR

Formalin-fixed, paraffin-embedded Rat Oviduct stained with Multi Cytokeratin Mouse Monoclonal Antibody (C11).

SDS-PAGE Analysis of Purified Multi Cytokeratin Mouse Monoclonal Antibody (C11).



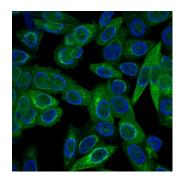

100 80 40 20 0 10¹ 10² 10³ 10⁴ 10⁵ FITC-A

Immunofluorescence Analysis of HeLa cells labeling Multi Cytokeratin with Multi Cytokeratin Mouse Monoclonal Antibody (C11) conjugated with CF640R(Green). The nuclear counterstain is DAPI (Blue)





Formalin-fixed, paraffin-embedded human Colon Carcinoma stained with Multi Cytokeratin Mouse Monoclonal Antibody (C11).

Specificity & Comments

This MAb recognizes cytokeratin 4, 5, 6, 8, 10, 13, and 18. This is a broad-spectrum antibody, which has been reported to differentiate epithelial tumors from non-epithelial tumors. Many studies have shown the usefulness of keratins as markers in cancer research and tumor diagnosis.

Known Applications & Suggested Dilutions

Flow Cytometry (1-2ug/million cells) | Immunofluorescence (1-4ug/ml) | Western Blot (1-2ug/ml) | Immunohistochemistry (Formalin-fixed) (1-2ug/ml for 30 min at RT)(Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95°C followed by cooling at RT for 20 minutes) | Optimal dilution for a specific application should be determined.

Limitations and Warranty

This antibody is available for research use only and is not approved for use in diagnosis.

There are no warranties, expressed or implied, which extend beyond this description. Company is not liable for any personal injury or economic loss resulting from this product.

Supplied As

200ug/ml of Ab Purified from Bioreactor Concentrate by Protein A/G. Prepared in 10mM PBS with 0.05% BSA & 0.05% azide. Also available WITHOUT BSA & azide at 1.0mg/ml.

Storage and Stability

Antibody with azide - store at 2 to 8° C. Antibody without azide - store at -20 to -80°C. Antibody is stable for 24 months. Non-hazardous. No MSDS required.

